



piccolo^{xpress}[®] chemistry analyzer

FAST, EASY, ACCURATE ON-SITE CHEMISTRY TESTING

REAL-TIME CLIA-WAIVED CHEMISTRY

EXCEPTIONALLY VERSATILE

Wide range of **CLIA waived** panels to increase your on-site diagnostic capabilities

- **Comprehensive Metabolic Panel**
ALB, ALP, ALT, AST, BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, TBIL, tCO₂, TP
 - **Basic Metabolic Panel**
BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂
 - **Lipid Panel**
CHOL, CHOL/HDL*, nHDLc*, HDL, LDL*, TRIG, VLDL*
 - **Lipid Panel Plus**
ALT, AST, CHOL, CHOL/HDL*, nHDLc*, GLU, HDL, LDL*, TRIG, VLDL*
 - **Liver Panel Plus**
ALB, ALP, ALT, AMY, AST, GGT, TBIL, TP
 - **General Chemistry 6**
ALT, AST, BUN, CRE, GGT, GLU
 - **General Chemistry 13**
ALB, ALP, ALT, AMY, AST, BUN, Ca, CRE, GGT, GLU, TBIL, TP, UA
 - **Electrolyte Panel** - Cl⁻, K⁺, Na⁺, tCO₂
 - **Kidney Check** - BUN, CRE
 - **Renal Function Panel**
ALB, BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, PHOS, tCO₂
 - **MetLyte 8 Panel**
BUN, CK, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂
-
- **Hepatic Function Panel**
ALB, ALP, ALT, AST, DBIL, TBIL, TP
 - **Basic Metabolic Panel Plus**
BUN, Ca, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂, Mg, Lactate Dehydrogenase
 - **MetLyte Plus CRP**
BUN, CK, Cl⁻, CRE, GLU, K⁺, Na⁺, tCO₂, CRP
 - **BioChemistry Panel Plus**
ALB, ALP, ALT, AMY, AST, BUN, Ca, CRE, CRP, GGT, GLU, TP, UA
 - **MetLac 12 Panel**
ALB, BUN, Ca, Cl⁻, CRE, GLU, LAC, K⁺, Mg, Na⁺, Phos, tCO₂

CLIA WAIVED

EXCEPTIONALLY FAST

From sample to complete panel results in 3 simple steps and approximately 12 minutes



1 / Add sample



2 / Insert disc



3 / Read results

EXCEPTIONALLY EASY

Fully automated system requires no special operating skills

- 0.1cc whole blood, serum or plasma
- Self calibrates with every test run - no manual data entry
- Full-color touch screen with intuitive user interface
- Results are printed on a self-adhesive hard copy or via electronic transfer, allowing immediate access to critical information at the point-of-care
- Large memory capacity stores up to 5000 patient records
- Externally accessible bi-directional USB data transfer technology facilitates EMR connectivity

EXCEPTIONALLY ACCURATE

Advanced technology ensures sound results

- Centrifugal action mixes the reagents and sample in a disposable self-contained disc – reducing procedural errors
- Accuracy and precision comparable to larger, more expensive laboratory analyzers
- iQC[®] system performs continuous quality control checks on every run, eliminating daily control testing requirements
- Every self-contained reagent disc is barcoded for automatic self-calibration

Distributed by:




50 Valleywood Drive
Markham, Ontario L3R 6E9
800.387.9643
800.268.1150 Français
905.470.2381 Fax
inter-medico.com



3240 Whipple Rd., Union City, CA 94587 | 800.822.2947 | www.piccoloxpress.com

©2015 Abaxis, Inc., Piccolo Xpress[®] and iQC[®] are registered trademarks of Abaxis, Inc. 888-3010 Rev.K

*Calculated

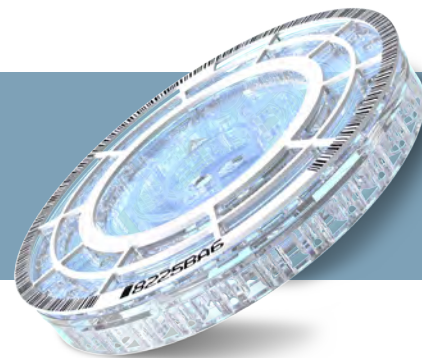

PICCOLO XPRESS PANEL DETAIL

Test	CLIA WAIVED												MODERATELY COMPLEX				
	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	
ALB	●		●		●		●						●			●	●
ALP	●		●		●								●			●	
ALT	●		●	●	●				●				●			●	
AMY			●		●											●	
AST	●		●	●	●				●				●			●	
BUN	●	●		●	●	●	●					●	●	●	●	●	●
Ca	●	●			●		●						●			●	●
CK												●			●		
Cl ⁻	●	●					●			●	●	●	●	●		●	●
CRE	●	●		●	●	●	●					●	●	●	●	●	●
GGT			●	●	●											●	
GLU	●	●		●	●		●		●			●	●	●	●	●	●
K ⁺	●	●					●			●	●	●	●	●		●	●
Na ⁺	●	●					●			●	●	●	●	●		●	●
PHOS							●									●	●
TBIL	●		●		●								●				
tCO ₂	●	●					●			●	●	●	●	●		●	●
TP	●		●		●								●			●	
UA					●											●	
CHOL								●	●								
CHOL/HDL*								●	●								
HDL								●	●								
LDL*								●	●								
nHDLc*								●	●								
TRIG								●	●								
VLDL*								●	●								
CRP														●	●		
DBIL													●				
LAC																	●
LDH												●					
Mg												●					●

* Calculated